

Fanttest

CE

Follicle-Stimulating Hormone (FSH) Rapid Test



Methodology

Using the principle of double antibody sandwich method.



Operation

Simple to operate and easy to interpret.



Detection

The detection is fast. The results can be interpreted in 3 minutes.
Detection level is no less than 25mIU/ml



Accuracy

High sensitivity and specificity



Hangzhou Fanttest Biotech Co., Ltd.

info@fanttest.com

+86 571 8633 7555

www.fanttest.com


Product background

The testing for FSH can help determine whether a woman is in the perimenopause stage. If a woman knows she is perimenopausal, she can take the appropriate steps to keep her body healthy and avoid the health risks associated with menopause, which include osteoporosis, increased blood pressure and cholesterol, and increased risk of heart disease.

Intended use

The Follicle-Stimulating Hormone (FSH) Rapid Test is a rapid chromatographic immunoassay for the qualitative detection of Follicle-stimulating hormone (FSH) level in urine to evaluate the onset of menopause in women.

Operation steps and result interpretation

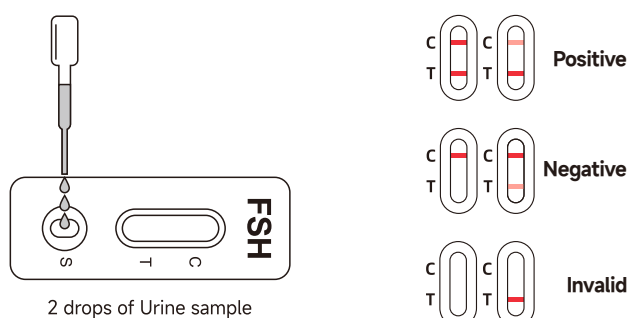

- **POSITIVE:** Two lines are visible and the test line (T) is the same as or darker than the control line (C). A positive result indicates that FSH levels are higher than normal and the subject you may be experiencing perimenopause.
- **NEGATIVE:** Two lines are visible, but the test line (T) is lighter than the control line (C), or there is no test line. A negative result indicates that the subject is probably not experiencing perimenopause in this cycle.
- **INVALID:** Control line (C) fails to appear. Insufficient specimen volume or incorrect test performance are the most likely reasons for an invalid result. Review the procedure and repeat with a new test. If the problem persists, discontinue using the test kit immediately and contact your local distributor.

Product information

Product name	Specimen	Format	Pack/box	Shelf life	Storage temperature	Certificate
Follicle-Stimulating Hormone (FSH) Rapid Test	Urine	Cassette	25T	24months	2-30°C	CE